Dr.Rock

Graphene Fast-Heating Red Light Waist Belt
石墨烯 速熱紅光腰帶

User Manual
使用說明書

ENGLISH 1 ~ 7

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Wish you health and peace.
願健康與平安歸給您。

⚠ Danger

Before using this device, please carefully read the manual and comply with the operating guidelines of the "Instructional Manual"; the company will not be responsible for the consequences of improper use or non-compliance with the given instructions for the correct use of particular product. This product is a personal care massage device and cannot be used as a substitute for medicines and medical devices.

⚠ Caution


- Do not use this device near to the water or moist use in the flammable.
- Do not use this device in the high electromagnetic environment, high temperature environment or in humid conditions.
- Do not use this device in the flammable and explosive places e.g. gas stations, high dust workshops etc.
- Do not disassemble the device by yourself. If the device failure occurs, please contact the after sales service to deal with the failure or to do its maintenance.
- Do not plug in or unplug the USB charging cable with the wet hands in order to avoid the danger of electric shock or any other kind of damage.
- Do not use this device with the other medical devices. Do not use more than 1 units of this product at the same time, otherwise it may cause the physical discomfort.
- Do not use it while driving a car, as it may cause physical discomfort due to excessive stimulation and that might lead to an accident.
- Do not use for any other purpose other than the message as it is clearly stated in the instruction manual to avoid accidents and any kind of malfunctioning.
- If any kind of discomfort is caused by the excessive stimulation during its use, please stop using the device immediately.
- This device should not be used by people those who are physically incompetent, have sensory or mental impairment, lack of relevant experience.
- Please do not give this device to the children or people without the ability to express consciousness to avoid accidents or physical discomfort.
- This device has heat, vibration and infrared red light function, please use with caution if you are not sensitive to the high temperature, vibration and infrared light intensity.
- Do not use the device in the following conditions: pregnant women, high fever, malignant tumors, untreated heart surgery, acute fibrosis, or other physical diseases.
- Do not use the device in the following conditions: program system, high fever, malignant tumors, untreated heart surgery, acute fibrosis, or other physical diseases.
- This device contains lithium battery, please do not disassemble, do not modify and do not throw into fire in order to avoid danger.
- When disposing of this device, please dispose of it in accordance with the relevant environmental protection regulations of the city.

Product Parameters

Device Name	Graphene Fast-Heating Red Light Waist Belt	Device Model	HXR-4.13D
Device Size	1152mm*222mm*14mm	Rated Power	60W
Battery Capacity	11.1V 2400mAh	Charging Voltage	5V, 2A
Device Weight	800g	Charging Port	USB Type-C

*Note: Because the belt has good elasticity, so there may be a 50% error between the measurement size and the actual size.

What is included in the package?

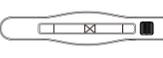Main Device*1

Type-C Charging Cable*1


Extension belt*1

Instructions Manual*1


Function and name of each part

Front side of the device

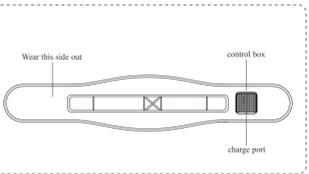Wear this side out control box charge port


Back side of the device

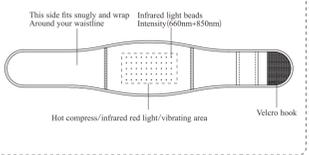This side fits snugly and wrap around your waistline Infrared light beads (Intensity:660nm+850nm) Hot compress/infrared red light/vibrating area Velcro hook

Device Operating Instructions

Control box panel

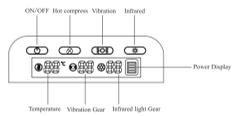
- Long press** button to turn on, the LED light before the button will be turned on into the standby mode, if there is no other operation within 2 minutes, the device will be Automatically turned off.
- During the Standby mode**, then short press the hot compress button or vibration button or infrared light button to operate the device, the default working time is 15 minutes, so after 15 minutes the device will be automatically turned off. If the user need to continue to use, then the user need to restart the device.
- Hot compress totally has 4 levels of heat, the temperature is: 40℃, 43℃, 45℃, 50℃, 53℃, 55℃, 58℃, 60℃. Short press the hot compress key to turn onto fiber heat gear.
- Vibration mode totally has 4 levels, the strongest mode, the relaxation mode, the working mode, the exercise recovery, by short pressing the vibration button key to turn onto fiber vibration mode. Infrared light totally has the 4 different levels, respectively: red light 1st gear one intensity level, red light 2nd gear two intensity level, red light 3rd gear three intensity level, red light 4th gear highest intensity level, short press the red light key to turn onto the intensity level of the infrared light.
- Long press** button to turn off the device.

Devic charging

If it is used for the first time or has not been used for a long time, please charge it as follows:

- The USB cable in the accessory should be plugged into charging port which is presently on the bottom left side of the device.
- Insert the USB cable into the power adapter and plug in the power adapter into the power socket.

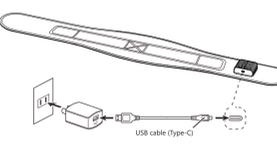- When the charging starts, the charging light effect appears on the power display of the device control box.
- The device will be fully charged after 4-5 hours of continuous charging, after fully charged battery the device light indicator will indicate that the battery is fully charged. Unplug the charging cable from the device, the device will be ready to use.

*Note: Do not use this device while charging.

Troubleshooting

Possible conditions	Possible causes/handling methods
Cannot start after pressing the power button	Please check whether the device has the sufficient power or not. Please confirm whether the device start-up operation is standardized (press and hold the key for more than 1.5 seconds). The device will automatically stop working after 15min of use, which is a normal phenomenon. If you need to continue to use the device then need to turn on the device again.
Stopped during its use	Please confirm whether the device is in low battery status, if the battery power is insufficient, please charge it and use it again after its fully charged. If the ambient temperature is particularly high or low, the device will not start, please charge at 10-15℃.
Unable to charge or very short life span	Please check whether the USB connection cable and the host power adapter are in good contact. The first charge or long time no use, please first full charge before its use. The device cannot be charged during stamp, please turn on the power before charging.
Heating/Vibration/Infrared function error	Please make sure the device is fully charged and switched on. Make sure that the corresponding mode of the control box is on one of A1-A2-A3-A4.
No sensation of hot compress function	Please check if your clothes are too heavy, it is recommended to use it after taking a shower out. Please make sure the heating function is on, and increase the heating level appropriately.
Rashness of the skin after use	This is a normal phenomenon after the skin has been cold heated, and will return to normal within half an hour after you stop using the device. Please check if your skin is allergic to the material of the device.
Tringling sensation	Please make sure there are no wounds on the skin surface, if there are wounds, please wait until the wounds have recovered before using the device.

*Note: If the problem still exists after following the above steps, please cut off the power and contact customer service to deal with it.

Maintenance and Storage Method after its use

● Maintenance method

- Be sure to unplug the power cable before its maintenance.
- Do not use detergent, thinner, gasoline or other chemicals to wipe or clean the device.
- When the device is dirty, it can be wiped with a dry soft cloth or dry towel to avoid washing in running water or other liquids.
- Please do not disassemble, repair or modify this device without the permission to avoid accidents or any kind of other possible failures.

● Storage Method

- Be sure to unplug the power supply during storage.
- Please do not keep it in a damp place or a place exposed to direct sunlight.
- If the device is not used for a long time, please fully charge the device and put it in the storage bag, then put it in a dry and ventilated place.
- Do not store it in a place where insects or pets can touch it.
- Please do not place it on cloth, paper products or wooden furniture that are easy to fade.

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